

Accessories, Interchangeable Blocks

For use with Techne® block heaters. Manufactured from anodised aluminium and all with a separate hole to accommodate a thermometer if desired. All blocks have dimensions (d x w x h) 95 x 76 x 51 mm and can be used in any combination (except 96-well blocks).

Aluminium Insert Blocks

Product Code	Tube Size (diam)	Number of Holes	Size d x w x h (mm)
F3501	Plain Block	None	95 x 76 x 51
F3502	6mm	30	95 x 76 x 51
F3503	10mm	20	95 x 76 x 51
F3504	12mm	20	95 x 76 x 51
F3505	13mm	20	95 x 76 x 51
F3506	15mm	12	95 x 76 x 51
F3507	16mm	12	95 x 76 x 51
F3508	19mm	8	95 x 76 x 51
F3509	25mm	6	95 x 76 x 51
F3510	10mm cuvettes	2 channels	95 x 76 x 51
F3512	Plain block	thermometer hole only	95 x 225 x 51
F4460	Plain block	thermometer hole only	95 x 76 x 51
F4461	7 and 9 mm	20/10	95 x 76 x 51
F4462	24mm	6	95 x 76 x 51
F4463	26mm	6	95 x 76 x 51
F4464	microcentrifuge 1.5ml tubes	20	95 x 76 x 51
F4465	0.5ml microtubes	30	95 x 76 x 51
F4466	Plastic spacer	None	95 x 37 x 51
F4467	Hi-Temp 96 block	96	95 x 151 x 61
F4468	Falcon round bottom plate block	96	95 x 151 x 61
F4469	Falcon flat bottom plate block	96	95 x 151 x 61
F4470	microcentrifuge 2.0ml tubes	20	95 x 76 x 51
F4471	0.2ml microtubes	72	95 x 76 x 51
F4473	Block for 96 x 0.2ml microtubes	96	95 x 151 x 61
F4474	1.0ml Porvair Plate	96	95 x 151 x 48
F4476	Block for Gelation Timer	1 sample cup	95 x 76 x 51



Aluminium Blocks


DBsoft Software DB-2D/3D/4D

DBsoft has been designed for use with Dri-Block® models DB-2D, DB-3D and DB-4D with the RS232 option. DBsoft is a Windows® based program that connects via an RS232 connection to your computer.

Programs can be created, saved and loaded at any time as well as running in real time mode. Definable parameters include time, temperature and ramp rates. Temperatures within the range of ambient to 200°C can be programmed.

DBsoft enables you to:

- Create, open, save and print programs
- Open, save, print and view logged data as a line graph
- Open, print and view logged data in a text format
- Create calibration certificates
- Run a program in real-time mode
- Send programs to the Dri-Block®
- Retrieve programs and test results from the Dri-Block®
- Erase test results
- Password protect programs

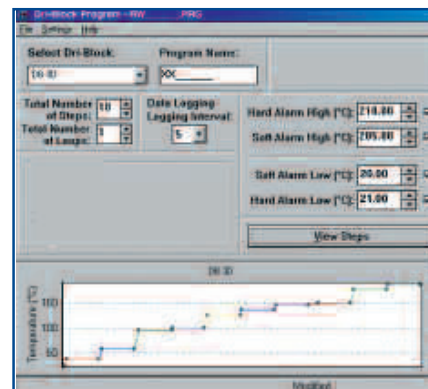

DBsoft

Ordering Information

Product Code	Description
FDBSOFT	DBsoft software for DB-2D/3D/4D