

## Features of TitroLine® 6000, 7000 and TITRONIC® 500

Features	TITRONIC® 500	TitroLine® 6000	TitroLine® 7000
pH/mV measuring input with reference input	–	1	1
Wireless electrode recognition	–	–	✓
Dead stop measuring input (2 x 4 mm socket)	–	✓	✓
Temperature measuring input (2 x 4 mm socket)	–	✓	✓
Interfaces	3 USB, 2 RS232	3 USB, 2 RS232	3 USB, 2 RS232
Titration to mV and pH endpoint	–	2 EP	2 EP
Dynamic and linear titration to equivalence point (EQ), mV and pH	–	1 EQ	2 EQ
Especially suitable for non aqueous titration	–	–	✓
Deadstop titration	–	✓	✓
pH-stat titration	–	–	✓
Manual titration	✓	✓	✓
Dosing tasks	✓	✓	✓
Preparation of solvents (manual or automatic with connected balance)	✓	–	✓
Formulae for miscellaneous calculations	✓	✓	✓
Standard methods	✓	✓	✓
Number of methods	5	> 10	50

## Ordering information

Set	Scope of delivery	Order no.
<b>T 500-M1</b> TITRONIC® 500 with magnetic stirrer	TITRONIC® 500 basic unit, magnetic stirrer TM 235, stand rod TZ 1510, electrode clamp Z 305, manual controller TZ 3880, power supply 100–240 V	285220210
<b>T 500-M2/20</b> TITRONIC® 500 with magnetic stirrer and 20 ml exchangeable unit	TITRONIC® 500 basic unit, magnetic stirrer TM 235, 20 ml exchangeable unit WA 20, stand rod TZ 1510, electrode clamp Z 305, manual controller TZ 3880, power supply 100–240 V	285220220
<b>TL 6000-M1/10</b> TitroLine® 6000 with 10 ml exchangeable unit	TitroLine® 6000 basic unit, magnetic stirrer TM 235, 10 ml exchangeable unit WA 10, incl. brown glass bottle for titration reagent, GL 45 and S 40 bottle adapters, hoses, drip off tube and titration tip	285220050
<b>TL 6000-M1/20</b> TitroLine® 6000 with 20 ml exchangeable unit	TitroLine® 6000 basic unit, magnetic stirrer TM 235, 20 ml exchangeable unit WA 20, incl. brown glass bottle for titration reagent, GL 45 and S 40-bottle adapters, hoses, drip off tube and titration tip	285220060
<b>TL 6000-M1/50</b> TitroLine® 6000 with 50 ml exchangeable unit	TitroLine® 6000 basic unit, magnetic stirrer TM 235, 50 ml exchangeable unit WA 50, incl. brown glass bottle for titration reagent, GL 45 and S 40 bottle adapters, hoses, drip off tube and titration tip	285220070
<b>TL 6000-M2/20</b> TitroLine® 6000 with 20 ml exchangeable unit and pH electrode	TitroLine® 6000 basic unit, magnetic stirrer TM 235, 20 ml exchangeable unit WA 20, incl. brown glass bottle for titration reagent, GL 45 and S 40 bottle adapters, hoses, drip off tube and titration tip. With pH-combination electrode IL-pHT-A120MF-DIN-N and buffer set.	285220080
<b>TL 7000-M1/10</b> TitroLine® 7000 with 10 ml exchangeable unit	TitroLine® 7000 basic unit, magnetic stirrer TM 235, 10 ml exchangeable unit WA 10, incl. brown glass bottle for titration reagent, GL 45 and S 40 bottle adapters, hoses, drip off tube and titration tip	285220140
<b>TL 7000-M1/20</b> TitroLine® 7000 with 20 ml exchangeable unit	TitroLine® 7000 basic unit, magnetic stirrer TM 235, 20 ml exchangeable unit WA 20, incl. brown glass bottle for titration reagent, GL 45 and S 40 bottle adapters, hoses, drip off tube and titration tip	285220150
<b>TL 7000-M1/50</b> TitroLine® 7000 with 50 ml exchangeable unit	TitroLine® 7000 basic unit, magnetic stirrer TM 235, 50 ml exchangeable unit WA 50, incl. brown glass bottle for titration reagent, GL 45 and S 40 bottle adapters, hoses, drip off tube and titration tip	285220160
<b>TL 7000-M2/20</b> TitroLine® 7000 with 20 ml exchangeable unit and pH electrode	TitroLine® 6000 basic unit, magnetic stirrer TM 235, 50 ml exchangeable unit WA 20, incl. brown glass bottle for titration reagent, GL 45 and S 40 bottle adapters, hoses, drip off tube and titration tip. With pH-combination electrode IL-pHT-A120MF-DIN-N and buffer set.	285220170

# TITRONIC® and TitroLine®

## New titrators and burettes



SI Analytix

[www.si-analytix.com](http://www.si-analytix.com)

SI Analytix

SI Analytix GmbH • Hattenbergstraße 10 • 55122 Mainz • Germany  
Telephone +49 (0)61 31-66 51 11 • Fax +49 (0)61 31-66 50 01  
[www.si-analytix.com](http://www.si-analytix.com) • [support@si-analytix.com](mailto:support@si-analytix.com)


# The new titration series

The new autotitrators TitroLine® 6000 and 7000 and the new piston burette TITRONIC® 500 make titration an easy task



## TitroLine® 6000 / TitroLine® 7000

Titrators

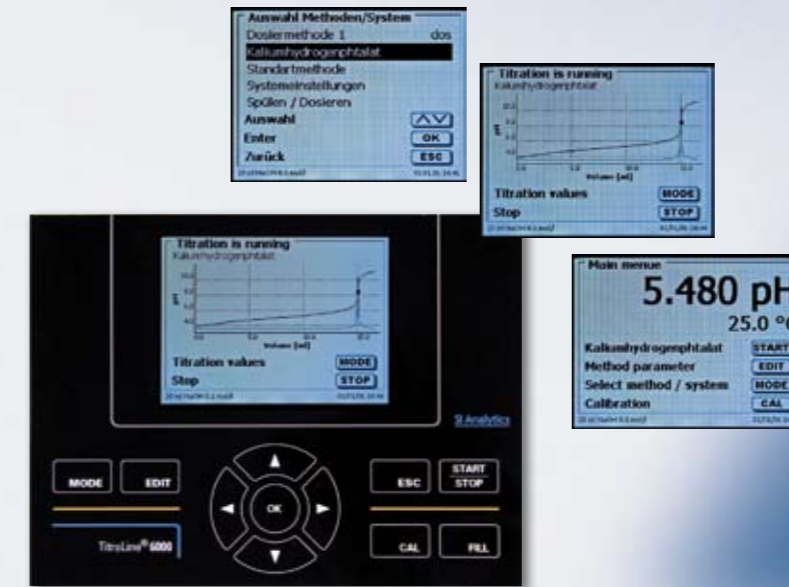

- ▶ The brilliant TFT display allows perfect readout even from the side.
- ▶ The intelligent exchangeable units are storing all relevant reagent and unit data.
- ▶ The automatic wireless recognition of SCHOTT® Instruments ID electrodes guarantees accurate calibration and measurements.
- ▶ The touch sensitive glass keypad is a pleasure to use.
- ▶ 3 USB and 2 RS-232 ports allow to connect e.g. USB keyboards, USB printers, USB storage media, stirrer, balance, PC and other SI Analytics instruments.



## TITRONIC® 500

Piston burette

